



CinemaVision

Setting New Standards
in Patient Relaxation




Making MRI Virtually Entertaining

Delivering unparalleled virtual reality, CinemaVision transports patients to a breathtaking world of immersive entertainment by combining fully-digital video and symphonic-quality audio. And, offering the ultimate in patient ease, CinemaVision offers a lightweight, adjustable integrated headset. Treat your patients to the highest quality multi-media experience and enjoy top-of-the-line performance!

CinemaVision Features Include:

- Widescreen video display with 16:9 aspect ratio or standard television 4:3
- Total resolution of 240,000 pixels on each display
- 45 degrees horizontal field-of-view
- Head mounting that produces the equivalent of a dramatic 62-inch screen viewed from 5.5 feet
- Proprietary MRI LaserLink technology that transmits audio-video signals from the console to the magnet room
- Digital audio system with 30dB noise-attenuating headset with 500Hz - 10KHz frequency response
- Compatible with 3T-rated field strength magnets, ensuring compatibility with all state-of-the-art MRIs

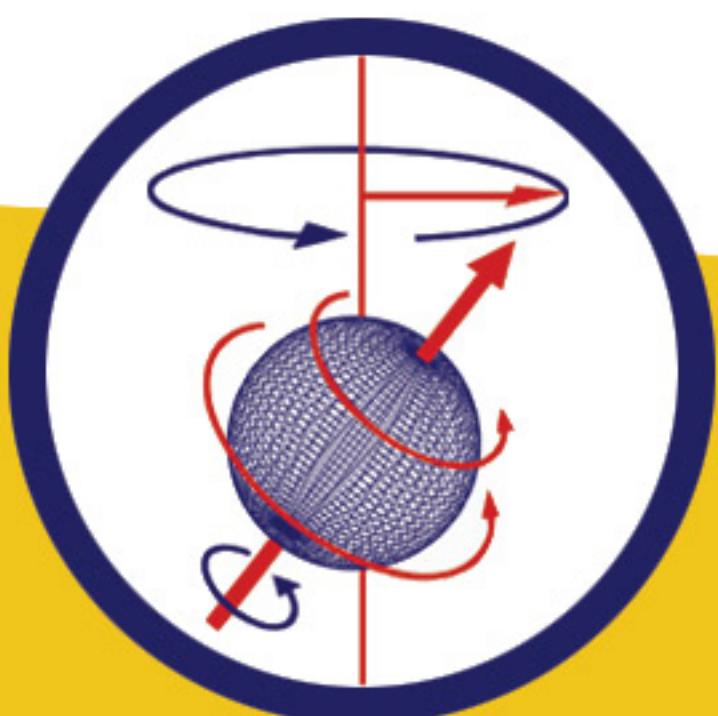

Resonance Technology, Inc.

818.882.1997 • 800.428.MRTV • www.mrivideo.com


CinemaVision

New Horizons in Patient Comfort

Transforming the patient's MRI experience, the CinemaVision is a versatile and immersive entertainment system. It provides multiple options from standard television, DVD, CD, AM/FM, and AUX inputs, while offering two-way communication and dramatically reducing MRI gradient noise. The head-mounted display fits completely within the head coil and operates seamlessly inside the magnet bore with no detrimental effects on signal/noise.



Eliminate appointment cancellations and sedation as patients enjoy a stunning sight and sound experience. With LaserLink technology, installation is just a matter of minutes. Patient set-up is quick and simple. Call today to learn more!

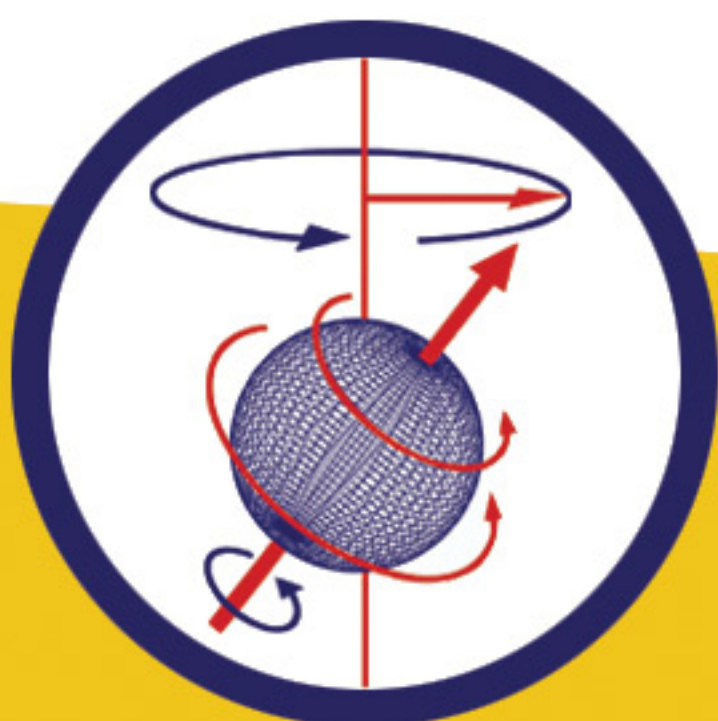

Resonance Technology, Inc.

818.882.1997 • 800.428.MRTV • www.mrvideo.com



Virtual Command Center ▲

The CinemaVision Controller manages entertainment choices for both the patient and technologist, enables two-way communication with patients, and more.

▲ Ultra-Sleek and Comfortable

CinemaVision's integrated audio-video headset, transducer and visor deliver top-quality entertainment and a virtual MRI comfort zone.

Handsome High-Tech Remote ▼

The CinemaVision remote control puts the functionality of the entire CinemaVision system at your fingertips!

